

EXPT.
NO 20

NAME

M T W T F S S


Page No.

YOUVA

Date

Aim: Determination of strength of Acid or Base by conductometer.

Principle: Solution of electrolytes conducts electricity due to the presence of ions. The specific conductance of a solution is proportional to the concentration of ions in it. The reaction between HCl and NaOH may be represented as,

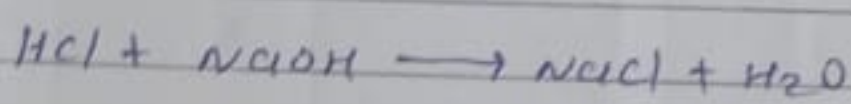

When a solution of hydrochloric acid is titrated with NaOH, the fast moving hydrogen ions are progressively replaced by slow moving sodium ions. As a result conductance of the solution decreases, this decrease in conductance will have place until the point is reached. Further addition of alkali raises the conductance sharply, as there is an excess of hydroxide ions. A graph is drawn between Volume of NaOH added and the conductance of solution. The exact end point is intersection of the two lines.

Teacher's Signature: _____

pH meter

Acid

Aim: To determine unknown concentration of acid by pH meter

Given data: Normality of NaOH = 0.2N volume of HCl = 10ml titration

$XN \text{ HCl} \rightarrow 0.2N \text{ NaOH}$

Pilot reading

micro reading.

Vme	pH	ΔpH
0	1.1	0.12
1	1.22	0.07
2	1.29	0.22
3	1.52	0.08
4	1.59	0.07
5	1.66	0.15
6	1.81	0.09
7	2.80 2.89	0.09
8	2.25	0.26
9	3.55	7.30
10	3.65	0.10
11	3.82	0.17
12	10.09	0.27
13	10.12 ₄	0.20
14	10.38	0.09
15	10.43	-

Vme	pH	ΔpH	ΔV	$\Delta pH / \Delta V$
0	2.15	0.05	1	0.05
1	2.2	0.15	1	0.15
2	2.35	0.13	1	0.13
3	2.48	0.13	1	0.13
4	2.66	0.32	1	0.32
5	2.98	0.27	1	0.27
6	3.25	0.43	1	0.43
7	3.68	1.24	1	1.24
7.2	4.52	0.06	0.2	0.3
7.4	4.68	0.14	0.2	0.7
7.6	5.12	0.23	0.2	1.15
7.8	5.35	0.13	0.2	0.65
8	5.48	0.20	0.2	1
8.2	5.68	2.34	0.2	11.7
8.4	8.02	0.19	0.2	0.95
8.6	8.21	0.14	0.2	0.7
8.8	8.35	0.10	0.2	0.5
9	8.45	0.07	1	0.07
10	8.52	0.03	1	0.03
11	8.55	0.13	1	0.13
12	8.65	0.07	1	0.07
13	8.75	0.40	1	0.40
14	9.15	0.13	1	0.13
15	9.3	-	-	-

...write
...what
...to do

EXPT. NO.	NAME	M T W T F S S						
		Page No.	YOUVA					
		Date:						

Result: The strength of the given strong acid is found to be $= 0.0823 \text{ N}$

The strength of the given strong base is found to be 0.08547 N

Teacher's Signature

[Signature]
16/10/23

EXPT. NO.	NAME	M T W T F S S	Page No.	YOUVA
			Date	

between a reference electrode and an electrode sensitive to the hydrogen ion activity when they are both immersed in the same aqueous solution. The reference electrode may be a Silver chloride electrode or a calomel electrode. The hydrogen-ion selective electrode is a standard hydrogen electrode.

Procedure :- $XN \text{ HCl} \rightarrow 0.1N \text{ NaOH}$
 $XN \text{ NaOH} \rightarrow 0.1N \text{ HCl}$

you have been given an unknown concⁿ on a 200ml of volumetric flask.

Dilute the given solⁿ to 200ml with distilled water.

Take some of this solⁿ in glass beaker and 20 ml of water into it.

Fill the burette with unknown concentration of HCl and NaOH.

Note down the reading by one ml interval.

Graph plot a graph between.

① $p^H \rightarrow V_{me}$ of titrant added

② $\Delta p^H / \Delta V \rightarrow V_m$ of titrant added


Teacher's Signature: _____

Sub/Ref.:


Date :
Time :

B.Tech CSE-E
IU2341230484

Scale: 1 cm = 1 unit

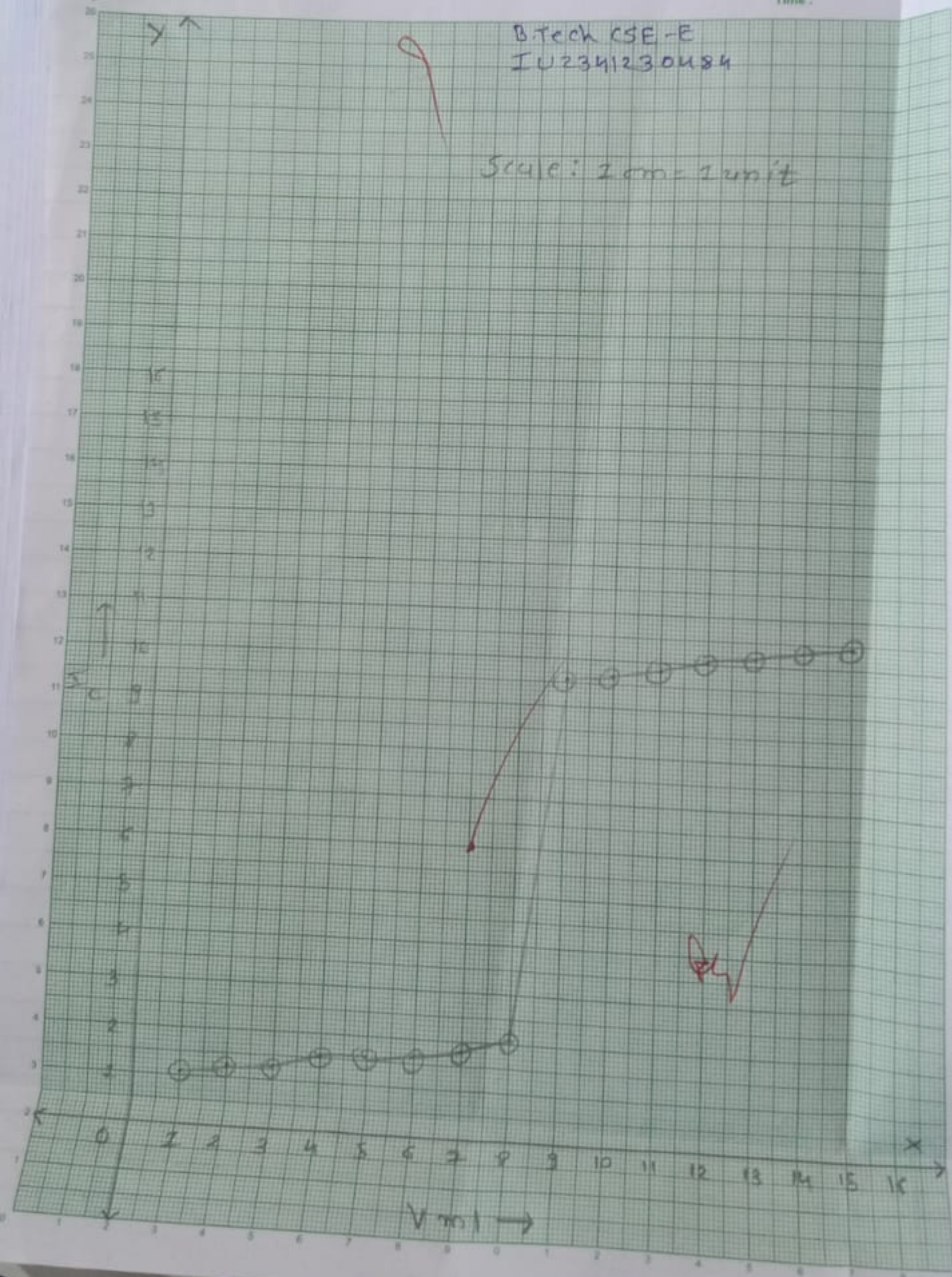

→ calculation:

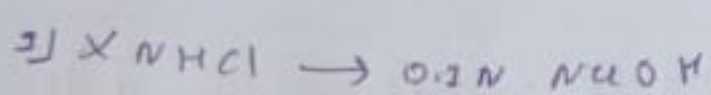


volume of sodium hydroxide $V_1 = \underline{10 ml}$

normality of sodium hydroxide $N_1 = \underline{0.2 N}$

volume of HCl $V_2 = \underline{8 ml}$

$$\begin{aligned} \text{normality of } HCl \ N_2 &= \frac{V_1 \times N_1}{V_2} \\ &= \frac{10 \times 0.2}{8} \\ &= 0.125 N \end{aligned}$$

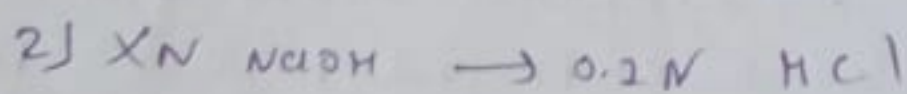

volume of ~~hydrochloric~~ hydrochloric acid $V_1 = \underline{10 ml}$

normality of ~~hydrochloric~~ HCl $N_1 = \underline{0.2 N}$

volume of $NaOH$ $V_2 = \underline{8}$

$$\begin{aligned} \text{normality of } NaOH \ N_2 &= \frac{V_1 \times N_1}{V_2} \\ &= \frac{10 \times 0.2}{8} \\ &= 0.125 N \end{aligned}$$